




Wi-Fi®


Wi-Fi set-up

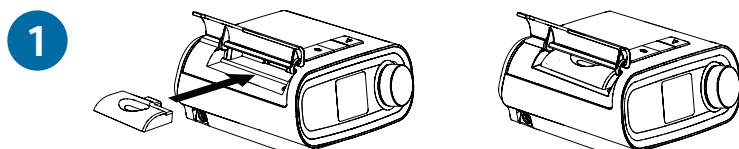

Verify that you have the DreamStation Wi-Fi Accessory Module

The Wi-Fi module is an optional accessory that may have been provided by the homecare provider.

- Open up the door on the left side of the device and identify if a module is inserted.
- If the module is gray in color then you have a Wi-Fi module attached to your PAP device.
- If the module is white then you have a Cellular accessory and no action is needed.
- If the compartment is empty then no action is needed.

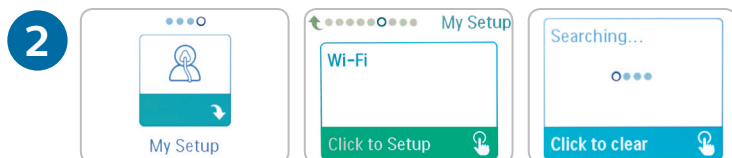

Finding your network

- In the PAP menu scroll to “My Setup” and press the control knob.
- In the “My Setup” menu scroll until you get to “Wi-Fi” and press the control knob.
- The PAP device will now search and identify all the networks within range.



Connecting to the network

- Choose the network you would like to connect to and press the control knob.
- Enter your Wi-Fi network passkey by scrolling through the menu and choosing the letters numbers / symbols that are necessary.
- Once the passkey has been entered scroll to the ✓ and press the control knob.
- The Wi-Fi accessory module will take a few seconds to connect to the network.


Complete

If additional assistance is needed please contact Philips Respironics at: 844-871-8067 or dreamstation.support@philips.com.



Bluetooth®

Bluetooth set-up

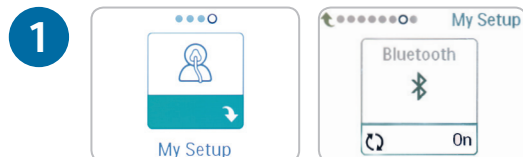

PAP

- In the menu screen scroll to “My Setup” and press the control knob.
- In “My Setup” scroll to “Bluetooth” and press the control knob.
- Confirm that your PAP device has Bluetooth turned “on.” Bluetooth will currently say “on.” If not, turn the control knob until the screen shows “Bluetooth On” and then press the button again to save the change.

2 Smartphone or tablet

- Enable the Bluetooth setting on your phone or tablet.

3 Smartphone or tablet

- The phone should now find the PAP Bluetooth.
- Select the “PR BT XXXX” and pair it to your device.
(Confirm that the characters “XXXX” match the last 4 characters in your DreamStation PAP serial number – see the label under your PAP device.)

4 Smartphone or tablet

- The mobile device will ask if you want to pair with the PAP device. Select “Yes.”



PAP

- Once the phone tries to pair with the PAP device the screen on the PAP will ask if you want to pair it with the device. Select “Yes.” Rotate the control dial to “Yes” then press the control dial to save the selection.



DreamStation has built in Bluetooth for connecting to external accessories such as DreamMapper* therapy self management system. DreamMapper can be downloaded for free to your mobile device or accessed via your computer at www.philips.com/mydreammapper



The Bluetooth(trademark symbol) word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and the use of such marks by Philips Respironics is under license. Other trademarks and trade names are of their respective owners.

*To see which therapy devices are compatible with DreamMapper, visit: www.dreammapper.com/compatible.

© 2016 Koninklijke Philips N.V. All rights reserved. Specifications are subject to change without notice. Trademarks are the property of Koninklijke Philips N.V. (Royal Philips) or their respective owners.

www.philips.com



Caution: US federal law restricts these devices to sale by or on the order of a physician.

RRDPGH EB 05/06/16 MCI 4107542 PN 1128895
1010 Murry Ridge Lane, Murrysville, PA 15668
800-345-6443 • 724-387-4000